

270 988

E-Z-ON™ MODIFIED VEST

TRANSPORT CHILDREN IN A PRONE OR SUPINE POSITION



FAA - 2004-17184-2

U.S. Patent No. 5,660,445

The E-Z-ON™ MODIFIED VEST (model number 101M2) is dynamically tested and meets FMVSS213. Shown with optional push button closures (D-Ring closures and standard). It features shoulder adjusters that enable the vest to "grow" with the child or adjust for seasonal clothing, and crotch straps that keep the lower hip strap of the vest positioned over the hip.

Transport children in a prone or supine position securely

Certain physical disabilities (such as spinal injury, body cast, long leg cast or hip spica) require a child to "lie down" when riding in a vehicle. The Modified Vest enables these children to be transported safely and securely.

Great for long trips...for all children

Now a child can nap comfortably in the car, and parents can have peace of mind knowing that their child is still secured to the vehicle's seat belt.

Uses your vehicle's existing seat belts

No Modification or locking clips are necessary. Just put the vest on the child (adjust it for a snug fit the first time you use it), lay the child on the seat, and secure the vehicle's existing seat belts after threading them through the sides of the vest and the leg belt. Place pillows as shown for comfort and protection in case of a collision. Complete, detailed instructions are included with each E-Z-ON™ MODIFIED VEST.

ORDERING INFORMATION

For most children ages 5-12 Model No. 101M2 (Double D-Ring Closure standard, Push-Button closure optional)

Installation and Usage Instructions for the Model 101M-2 Modified E-Z-ON™ VEST

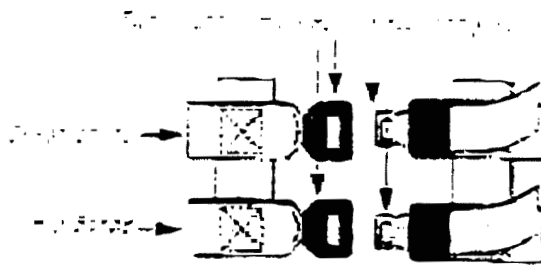


Diagram A

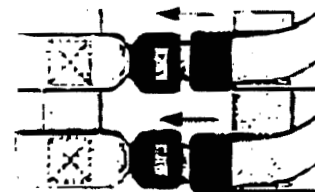


Diagram B

STEP #1: PUT THE VEST ON THE CHILD....

Place the Modified E-Z-ON™ VEST on the child with the closures in the front and the parachute emblem over the left shoulder strap.

If your vest has push-button closures (diagrams A & B):

- Insert the metal tongue into the push-button buckle on the chest strap until you hear a snapping sound. Tug to make sure the connection is secure.
- Repeat the above procedure on the hip strap.
- To adjust chest & hip straps:
 - To tighten: pull on the free ends of the belts.
 - To loosen: hold the adjusters at a right angle to the straps and pull.
- To unfasten chest & hip straps: press the button in the buckle.

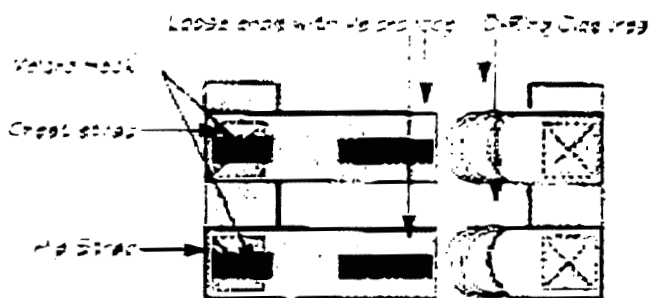


Diagram C

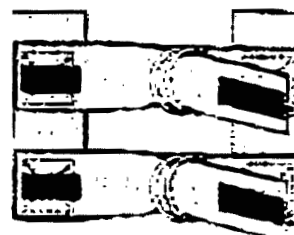


Diagram D

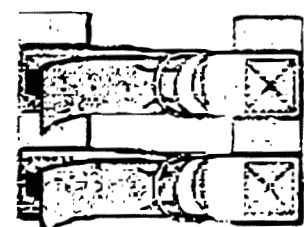


Diagram E

If your vest has D-Ring closures (diagrams C, D, E):

- Locate the loose end (with Velcro loop) of the chest strap (diagram C).
- Thread the loose end through both D-Rings on the opposite strap (diagram D), and pull the strap tight for a snug fit on the child.
- Then, loop the loose end over the first D-Ring and through the second D-Ring. Press Velcro loop on the loose end onto Velcro hook on the chest strap. (diagram E).
- Repeat the above procedure with the loose end of the hip strap.



Diagram E

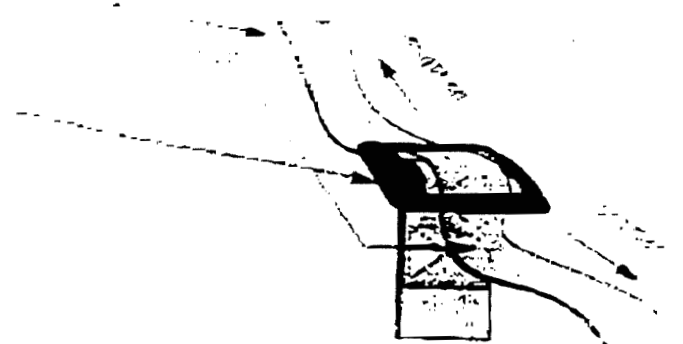


Diagram F

STEP #2: ADJUST THE FIT...

The shoulder straps should be adjusted so that the hip strap of the vest (the bottom strap with the push-buckle or D-Ring closure) is positioned over the pelvic bone when the vest is worn by the child.

To adjust the shoulder straps:

a. Pull the loose end of the shoulder strap (see diagram E) out of the bar slide. The bar slide and the loose end of the shoulder strap should now look like Diagram F.

b. To shorten the shoulder strap, hold the bar slide with one hand while pulling the loose end of the shoulder strap (see Diagram F).

c. To lengthen the shoulder strap, hold the bar slide with one hand while pulling the shoulder strap (see Diagram F).

d. Secure the shoulder strap by threading the loose end back through the bar slide (Diagram G). Pull loose end to tighten (Diagram H).

e. Repeat the above steps on the other shoulder strap.

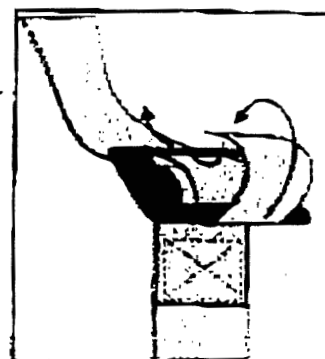


Diagram G

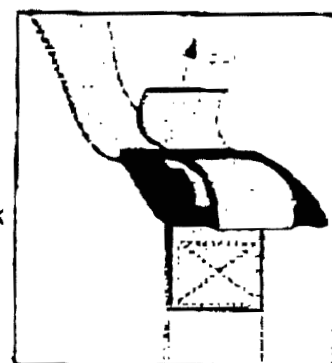


Diagram H

The crotch straps help keep the Modified Vest positioned correctly over the pelvic bone during use. To adjust the crotch straps:

a. With the vest on the child, secure the crotch strap fastener.

b. To shorten the crotch strap, hold the fastener at a right angle to the crotch strap and pull the loose end of the crotch strap (See Diagram I).

c. To lengthen, pull on the crotch strap (See Diagram I).

d. Repeat on other side.

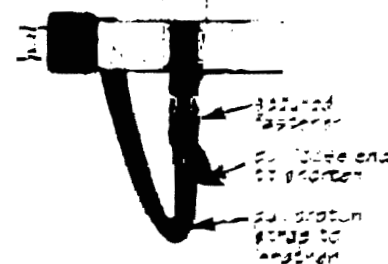


Diagram I

STEP #3: IMPORTANT!

PROTECT THE HEAD AND NECK

Place a small pillow, thick foam, rolled towel or rolled blanket between the seat back and the child's head. This will minimize neck flexion in the event of a collision.

BE SURE TO KEEP THE CHILD'S HEAD AND FEET AWAY FROM DOORS OR WALLS OF VEHICLE. Place a pillow between the door or wall and the child's head to help prevent possible injury during a lateral collision. The same precaution should be used at the child's feet

Please call
E-Z-ON PRODUCTS, INC. OF FLORIDA
if you...

- have questions;
- need assistance with installation;
- need additional information;
- or, want to place an order.

1-800-323-6598



E-Z-ON PRODUCTS, INC. OF FLORIDA

605 Commerce Way West • Jupiter, FL 33458 • Phone (561) 747-8920 • Fax (561) 747-8779

DEATH OR SERIOUS INJURY CAN OCCUR IF
MANUFACTURER'S INSTRUCTIONS ARE NOT
FOLLOWED CORRECTLY

E-Z (H) MOBILE HOME MOBILE
Use proper, proper, or lateral position
Mount (H) MOBILE HOME MOBILE when used with (2)
vehicle seatbelts

Age 2, 20-44 lbs
20" or less use one vehicle belt only through harness

Reduce slack in webbing

Mfg by E-Z (H) PRODUCTS, INC. OF FLORIDA

Jupiter, Florida (800) 373-6086 U.S. Pat. No. 5,000,445



DEATH OR SERIOUS INJURY CAN OCCUR IF
MANUFACTURER'S INSTRUCTIONS ARE NOT
FOLLOWED CORRECTLY

E-Z ON MOBILE HOME MOBILE 101642

Use proper, proper, or lateral position.

Mount (H) MOBILE HOME MOBILE when used with (2)
vehicle seatbelts

Age 2, 20-44 lbs.

20" or less use one vehicle belt only through harness

Reduce slack in webbing

Mfg by E-Z (H) PRODUCTS, INC. OF FLORIDA

Jupiter, Florida (800) 373-6086 U.S. Pat. No. 5,000,445

